JOEY Hardwire

BRUKSANVISNING (S)

USER MANUAL (GB)

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MODO DE EMPLEO (E)

MANUALE D'USO (I)

BRUKSANVISNING (N)

S

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The contents of these instructions must be read and understood before the product is taken into use. Keep the instructions easily accessible for further reference. Installation and repair must for security reasons and for warranty conditions be carried out by Autoadapt AB or by an Autoadapt AB authorized workshop. Joey wheelchair lift is safety tested. If the product is altered, or not installed according to the installation instructions, Autoadapt AB denies any liability.

1. PRESENTATION

The Joey lift offers an easy way of loading and unloading a wheelchair in a vehicle. The wheelchair is driven onto a platform where it is secured by three belts. To load the wheelchair the user only needs to press a button on the hand control. Joey is designed for vehicles with enough luggage area to fit the lift platform and rail, such as MPV and mini buses.

If you have a Joey platform lift installed in your vehicle make sure to inform your insurance company of the fact.

When Joey is to be removed from your vehicle this will be at an estimated cost of 1-3 hours, depending on the vehicle model and on how the installation has been made.

2. GENERAL SAFETY REGULATIONS

-  **Park the vehicle on horizontal, level surface**
-  **Always turn the I/O switch of the handcontrol OFF when Joey is not in use**
-  **If the Joey lift is not in use for more than a week the circuit breaker, placed on the right hand side of the cover, shall be switched to OFF. This prevents the battery from discharging.**
-  **Make sure no persons or objects are within the loading area when the lift is operated!**
-  **Keep hands and fingers away from lift tower during lift operation! Risk for pinch accidents**
-  **The wheelchair / scooter shall be secured to the platform according to instructions.**

3. OPTION


Retractable platform belts can be ordered as option. These belts are easier to use than the standard belts; pull out the hooks, mount them on the wheelchair docking plate and tighten by turning the handle clock wise.

Please note that the retractable belts demand more room between platform and vehicle rear door. The compartment space must be taken into consideration.

4. LOAD

1. Make sure the circuit breaker, to the right on the cover, is switched to ON.


(If the lift is not used for more than a week the circuit breaker shall be turned OFF to avoid a discharged battery).

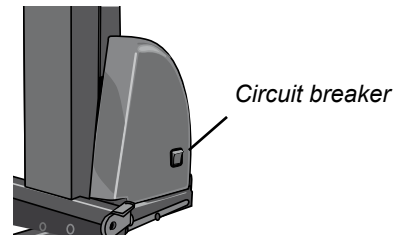2. Press the ON/OFF switch of the hand control to **ON**.


Make sure there are no people or objects within the loading area!

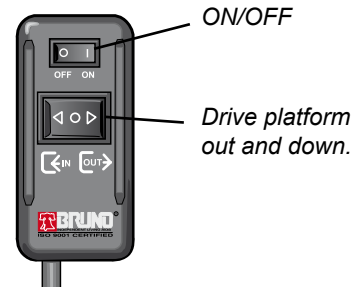
3. Drive the platform out and down by pressing the switch to **OUT**. Keep the switch pressed down until the platform has properly reached the ground.

4. Clear the platform and position the mobility device on the platform. Insert the hooks of the belts into the appropriate holes in the docking plate. Pull the straps to tighten! (see illustration)

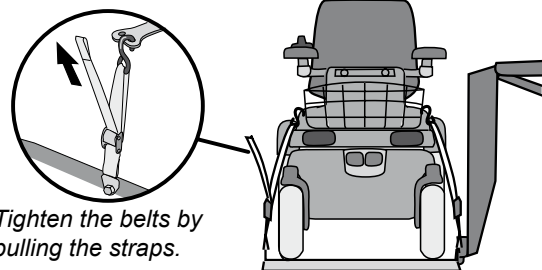5. Press the hand control switch marked IN, until the platform has reached end position inside the vehicle. Let go of switch when the lift has stopped.

The green LED on the lift vertical mast shines when the lift is in the inside, docked position.

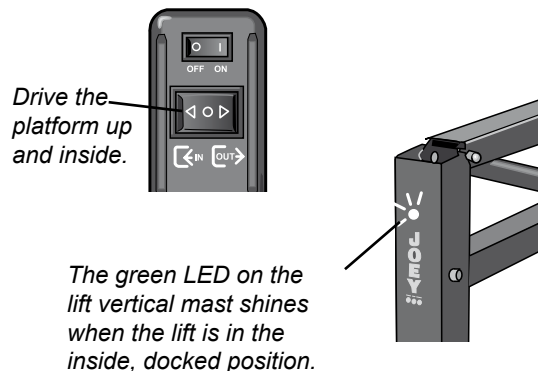
6. Press the hand control ON/OFF switch to **OFF** and place it in its holder on the wall.

5. UNLOAD

1. Press the hand control ON/OFF-switch to **ON**.


Make sure there are no people or objects within the loading area!

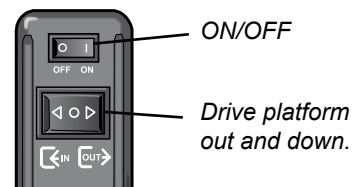2. Drive the platform out and down by pressing the switch to **OUT**. Keep the switch pressed down until the platform has properly reached the ground.

3. Release the securing belts by first pressing the strap release latch. Then dismount the hooks from the wheelchair docking plate.

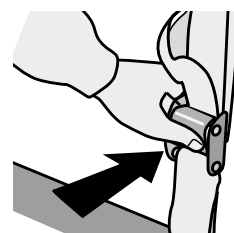4. Move the mobility device away from the platform.


5. Press the hand control switch marked IN, until the platform has reached end position inside the vehicle. Let go of switch when the lift has stopped.

The green LED on the lift vertical mast shines when the lift is in the inside, docked position.

Drive the platform up and inside

The green LED on the lift vertical mast shines when the lift is in the inside, docked position.

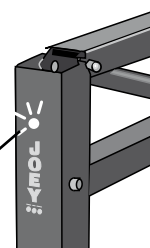6. Press the hand control ON/OFF switch to **OFF** and place it in its holder on the wall.

If the Joey lift is not in use for more than a week the circuit breaker on the cover shall be switched to OFF!

(To prevent the battery from discharging)


Circuit breaker

6. PINCHING POINTS

When the lift is being operated it is important to keep fingers away from pinching points (marked with a circle in the illustration).

Pinching points

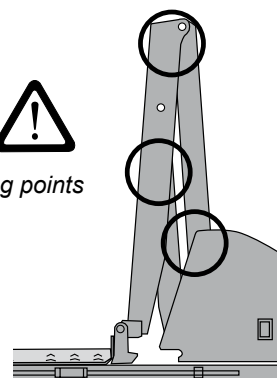7. EMERGENCY OPERATION

In cases of emergency the Joey lift can be manually operated. There are two actuators; one for IN/OUT and one for UP/DOWN.


1. Switch the circuit breaker to **OFF**!

2. Disengage the emergency tool from the two clips.

3. For moving up or down: insert the tool in the top entrance, according to illustration. Mount the tool to a drilling machine and rotate clockwise for UP, or counter clockwise for DOWN.


4. For moving in or out: mount the tool in the entrance facing the seats, according to illustration.

5. Fix the tool back in its place at the side of the lift.

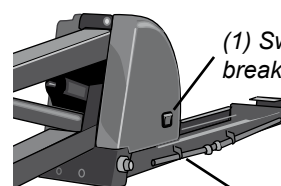(1) Switch the circuit breaker to OFF

(2) Emergency tool

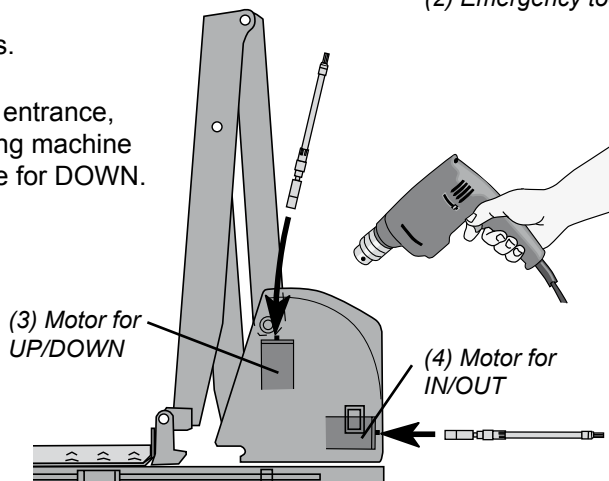(3) Motor for UP/DOWN

(4) Motor for IN/OUT

8. TROUBLE SHOOTING

PROBLEM	POSSIBLE CAUSE / SOLUTION
Lift does not operate.	<i>Verify that the hand control ON/OFF switch is ON.</i>
Lift operates slowly, lacks power or stalls.	<i>Weak vehicle battery. Charge battery. Load is too heavy. Do not exceed maximum load of 158 kg.</i>
The platform can be lifted up but not powered in.	<i>Adjustment of IN enable limit switch is needed for full horizontal retraction. Contact authorized dealer /technician.</i>
Lift makes a "machine gun" noise when platform is at lowest position.	<i>This is normal. The lift's actuator includes a clutch which has run up or down as far as possible.</i>

9. TECHNICAL SPECIFICATIONS

<p>Total length: 115 cm*</p> <p>Total width: 112 cm</p> <p>Total height: 97 cm</p> <p>Weight: 75 kg</p> <p>Platform measurements: 69 cm x 112 cm</p> <p>Maximum lift capacity: 158 kg</p> <p>12-Volt DC motor: draws less than 10 Ampere with load under 158 kg.</p> <p>* With standard securement belts.</p>	<p><i>The lift fulfils requirements according to EMC 95/54/EC, 2004/108/EC and is crash tested according to FMVSS 208S (in Europe ECE-R-17)</i></p> <p><i>CE-marked.</i></p> <p><i>The lift is equipped with a sensor that stops the lift if obstructed.</i></p> <p><i>The lift can be manually operated in cases of emergency.</i></p>
--	---

10. MAINTENANCE

Periodic maintenance and inspection

- Make sure that the green LED of the vertical mast shines as the lift has reached its inner (parked) position.
- Examine all electrical connections and remove corrosion.
- Check fuse connections by the vehicle battery and keep them clean.
- Visually examine lift; look for signs of structural damage (deformations or crackings). Contact authorized dealer for necessary reparations.
- Check that sensor for obstructions is functioning when lift is traveling in.
- Examine vehicle tires, suspension and response while driving.
- Review correct operation with user.
- Verify that pivot point fasteners are tight.

Annual maintenance and inspection

- Take vehicle with lift to an authorized dealer for inspection.
- Check that platform is horizontal (with load), x- and y-wise, both when travelling horizontally and when lifted up or down. Adjust if necessary.
- Grease pivot points (four places) with teflon based grease.
- Grease gear rack with multi-purpose grease.
- Test Joey battery.

Annual inspection record:

	Date	Mileage	Workshop/dealer	Comments
Service no. 1				
Service no. 2				
Service no. 3				
Service no. 4				
Service no. 5				
Service no. 6				
Service no. 7				

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